


Protocol for IHC-Fr(Anti Iba, Rabbit)

Wistar Rat (7 week old, ♂) was perfusion-fixed with 4% paraformaldehyde, replaced sucrose, and prepared 50µm brain section by microtome.

Wash : 0.3 % TritonX-100 in PBS、5 min×3
↓
Blocking : 1 % BSA, 0.3 % TritonX-100/PBS, 2hour , RT
↓
Primary antibody : rabbit anti-Iba1(1 : 1000), 1 % BSA, and 0.3 % TritonX-100 in PBS、
Over night 、4°C
↓
Wash : 0.3 % TritonX-100 in PBS、5 min×3
↓
Secondary antibody : AlexaFluor488 anti-rabbit IgG(1:1000), 1 % BSA, and 0.3 %
TritonX-100 in PBS, 2hour, RT
↓
Wash : 0.3 % TritonX-100 in PBS、5 min×3
↓
Mount

IHC data sample

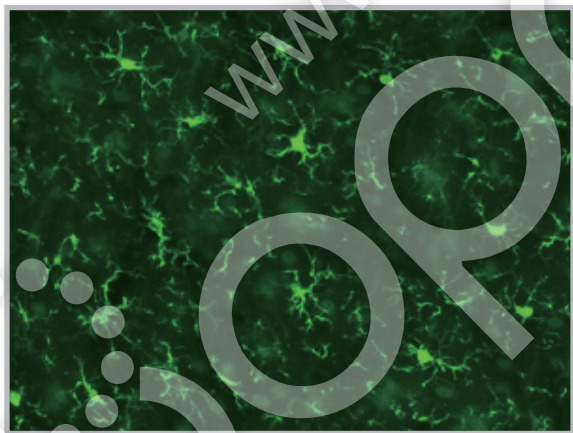

产品列表

产品编号	产品名称	包装
019-19741	小胶质细胞/巨噬细胞特异性蛋白抗体（免疫组化） Anti Iba1, Rabbit (for Immunocytochemistry)	50 µg

Wako

和光純薬工業株式会社
Wako Pure Chemical Industries, Ltd.

boppard
宝柏·中国

北京 Tel: 010 85804838
上海 Tel: 021 62884751
广州 Tel: 020 87326381
香港 Tel: 852 27999019
Email: info@boppard.cn

