

# **MONOCLONAL ANTIBODY**

For research use only, Not for diagnostic use.

## Catalog No. CAC-SHI-EXO-M02-EX

# Anti CD63 [ Clone : 8A12]

## **BACKGROUND**

CD63 (also known as LAMP-3, Melanoma-associated antigen ME491, TSPAN30, MLA1 and OMA81H) is a cell surface glycoprotein which belongs to the tetraspanin superfamily. CD63 is known to complex with integrins. CD63 is expressed on activated platelets, monocytes and macrophages, and is weakly expressed on granulocytes, T cell and B cells.

Product type Primary antibody

Immunogen Human CD63(a.a 104-202)

Rased in Mouse

Myeloma -

Clone number 8A12 lsotype lgG2a

Source Culture supernatant

**Purification** Affinity purified by Protein A

**Buffer** Phosphate buffered saline containing 0.1% NaN<sub>3</sub> as a preservative

Concentration1 mg / mLVolume100 ulLabelUnlabeledSpecificityHuman CD63

Cross reactivity Human. Other species is not tested.

Storage Store below - 40°C

Aliquot to avoid cycles of freeze/thaw..

Application notes

Western blotting: 1/500 - 1/2,000 (0.5 - 2ug/mL, non-reductive condition)

Recommended

**Immunoprecipitation**: 1/200 – 1/1,000 (1-5ug/mL)

dilutions


\* Experiment example : Serum and Culture supernatant Exosome

Other applications have not been tested.

Optimal dilutions/concentrations should be determined by the end user.

References -

#### **ANTIBODY CHARACTERIZATION**

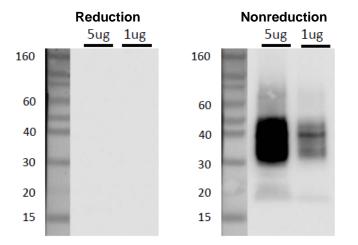

Sample:
 C32 cell Exosome (5ug, 1ug)
Primary antibody:
 anti-CD63 (8A12) (1ug/mL)
Secondary antibody:
 Anti-mouse IgG-HRP

Fig.1 Western blot analysis of CD63 in C32 cell Exosome with anti-CD63 (8A12) antibody



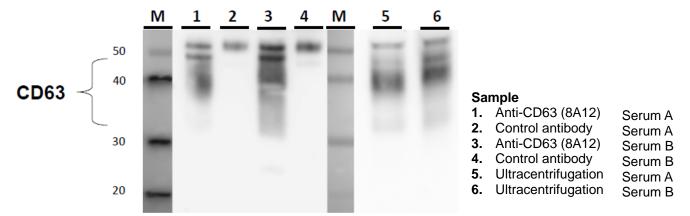


Fig.2 Immunoprecipitation in Human serum exosome

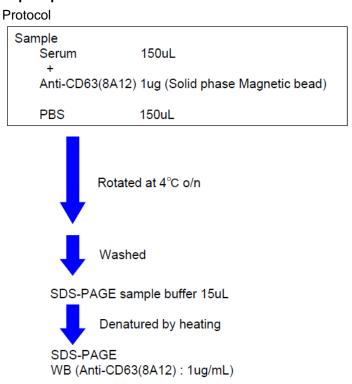


# www.cosmobio.co.jp

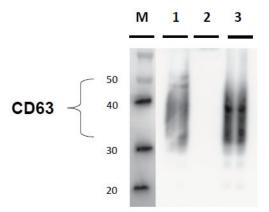

#### Sample

- 1. Anti-CD63 (8A12)
- 2. Control antibody
- C32 cell Exosome (0.4 ug)

Fig.3 Immunoprecipitation in C32 cell (Melanoma cell line) Exosome

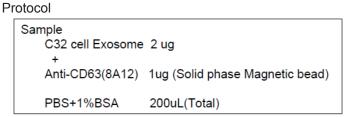



SDS-PAGE sample buffer 15uL



Denatured by heating

SDS-PAGE (3uL)

WB (Anti-CD63(8A12): 1ug/mL)

# **RELATED PRODUCT:**

Product Name	Application	Quantity	Cat#
Anti CD9 (12A12) Monoclonal Antibody	WB / IP / ELISA	100 ug	CAC-SHI-EXO-M01-EX
Anti CD63 (8A12) Monoclonal Antibody	WB / IP	100 ug	CAC-SHI-EXO-M02-EX
Anti CD81 (12C4) Monoclonal Antibody	WB / IP	100 ug	CAC-SHI-EXO-M03-EX



宝柏•中国

www. boppard.cn info@boppard.cn

北京 Tel: 010 85804838 上海 Tel: 021 62884751 广州 Tel: 020 87326381 香港 Tel: 852 27999019

For research use only, Not for diagnostic use.

Manufactured by SHIONOGI & CO., LTD.



COSMO BIO CO., LTD.

Inspiration for Life Science

TOYO 2CHOME, KOTO-KU, TOKYO, 135-0016, JAPAN

URL: http://www.cosmobio.co.jp e-mail: export@cosmobio.co.jp

[Outside Japan] Phone: +81-3-5632-9617 [国内連絡先] Phone: +81-3-5632-9610 FAX: +81-3-5632-9618 FAX: +81-3-5632-9619